


AFZAN LAB MANUAL

DIAGRAM

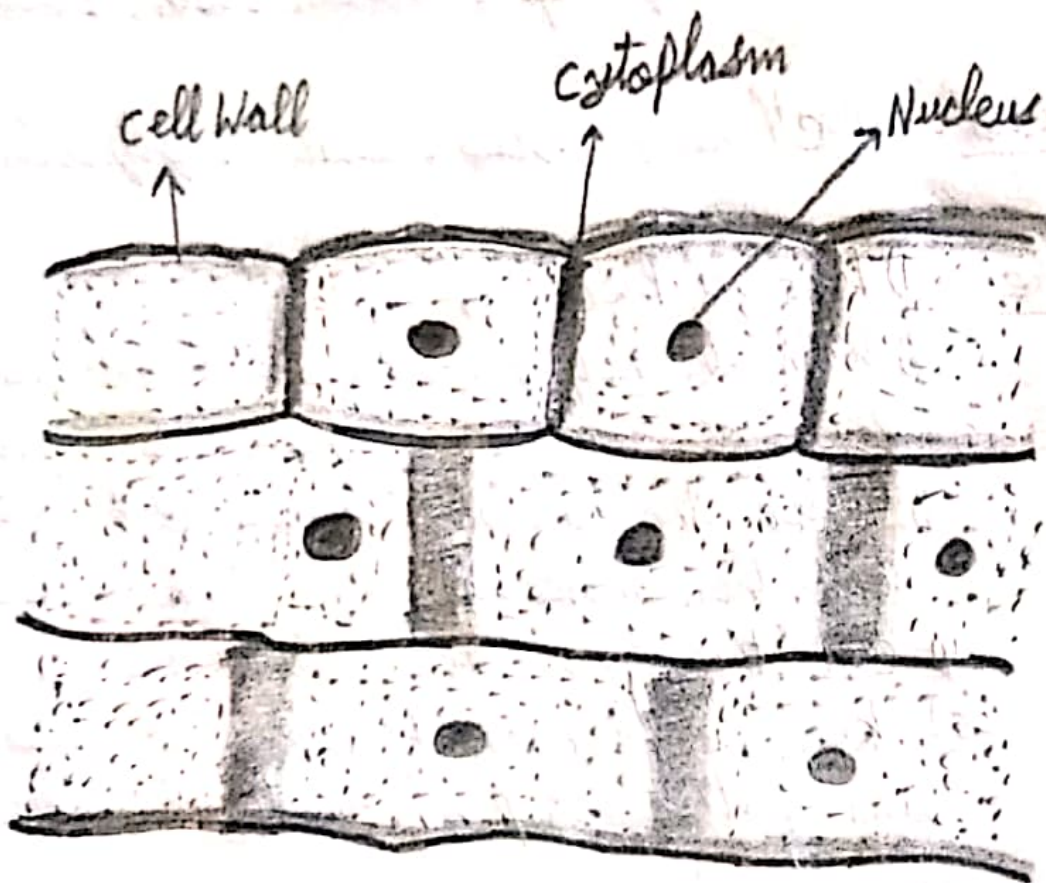

AFZAN LAB MANUAL

Experiment No. 2(A)

Page No.

Date 06/08/26

- To Prepare and study Plant cells from Onion Peel.
- Requirement - Microscope, cover slips, Watch glass, slide, blade, onion, forceps, Muslin cloth, blotting paper
- Chemicals - Methylene blue, water, glycerine.
- Procedure -
 - 1) Take a piece of Onion.
 - 2) With the help of forceps remove the white layer or Onion Peel from its concave side.
 - 3) Place the removed onion Peel into water and then make it dry with muslin cloth.
 - 4) Cut this Onion Peel into pieces.
 - 5) keep a piece of peel on the glass slide and put a drop of methylene blue stain on it.
 - 6) Put a drop of glycerine on stained peel and place a cover slip on it.
 - 7) Clean the slide from edges.
 - 8) Observe this mount first under low power and then under high power microscope.

Observations	Result.	
	Low Power Microscope	High Power Microscope
Observe the shape of the cell	- Shape of the cell will be rectangle	- Shape of the cell will be rectangular with distinct cell wall.

FZAN LAB MANUAL

Experiment No. 2(a)

Page No.

Date

- Observe the nucleus of the cell.	There is a deep colored central round body, which is nucleus.	Nucleolus is visible within the Nucleus.
- Observe the cytoplasm	-----	It is deep coloured on inner surface of cell wall.
- Observe the vacuoles	-----	Protoplasm contains three vacuoles which has its outer covering.

- Result - Onion Peel is made up of many cells having cell wall, nucleus, ~~vacuole~~ vacuoles and cytoplasm.